

## In Stock


## Only \$ 18,975

Usually \$26,500

Code: SANYMDF-DU702VX-PE



Offering 729 litres of storage (576 x 2" cryoboxes) and enhanced safety due to the dual cooling system.

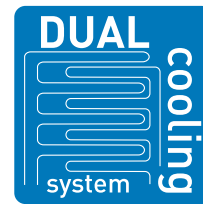

## As safe as you can get!

### Dual Cooling System

Two independent refrigeration systems provide the best redundancy you can have! – if one refrigeration circuit fails – the other circuit will maintain -70°C.

- Quiet in operation – only 52dB(A)
- Filter-less condenser design: eliminates the need for routine filter cleaning
- Suitable for operation on 220/240V 50Hz 10 amp supply
- Automatic vacuum relief port for quicker access to samples after door closure
- Hot gas bypass combined with fivefold door gasket assists in minimising frosting and ice build up around door seal
- Alarms for high and low temperature, door open and power failure
- Intelligent ECO Mode: optimising compressor operation to match the freezer's running condition and save electricity costs
- User friendly Colour Graphical Control Panel with digital display
- USB port for downloading logged data
- **VIP Plus insulation:** provides typically 24 % more storage capacity compared to conventionally insulated freezers with the same footprint
- Castors for mobility and locating feet when in working position
- Keyed door lock for security
- Dry remote alarm contacts
- Three access ports
- 5-year warranty