


## Next Generation Incubators for Optimum Cell Culture

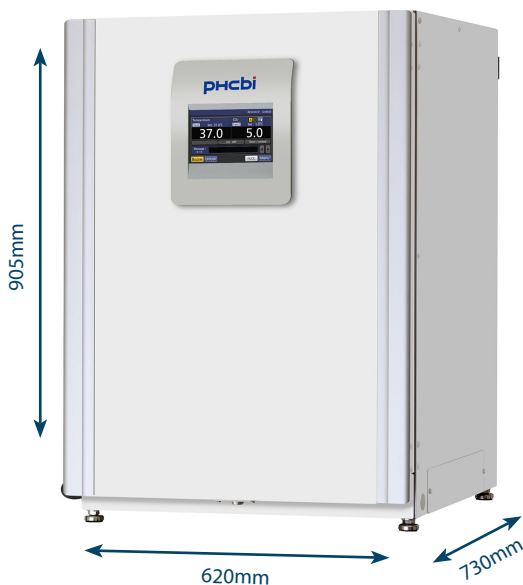

**Only \$7,555**

**Usually \$10,796**

Includes 4 shelves

Code: SANYMCO-170AICUVL-PE

*Price valid for current stock only*

**2 Units In Stock**

Superior culturing and with the ease of reactive contamination control via UV irradiation of the airstream for a programmable time period every time the door is closed and pro-active contamination control via the inCu-saFe™ corrosion-resistant blended copper alloy stainless steel interior!

- Inner glass door for viewing cultures
- Temperature: Ambient +5°C to 50°C (±0.1°C)
- CO<sub>2</sub>: 0 to 20% (±0.15%)
- Humidity: 95%RH (±5%RH)
- Direct heat and air jacket for environment stability
- Dual wavelength single beam dual detector IR CO<sub>2</sub> sensor with continuous auto zero calibration
- Programmable SafeCell™ UV system (isolated from cell cultures) irradiates airflow and humidity pan to prevent contamination without affecting cell cultures in-vitro
- Field reversible door
- Peltier powered dew stick to draw condensation away from inner glass door
- PID control system with colour LCD display
- USB port for data transfer
- Fully rounded corners and integrated shelf supports reduce contamination and simplify cleaning
- Integrated side wall tray supports – no bolts or screws that can harbour contamination
- 5-year warranty

Ref. Code: SANY-CO<sub>2</sub>

Valid Until: 31 Mar 2025

